
Air-conditioning Monitoring System

TAICOM-RP Monitor

Features of TAICOM-RP Monitor

- Continuous monitoring of manufacturing environment
- Support of scientific research for nonconformity and deviation management
- GAMP5 (including Part 11) compliant

Continuous monitoring of manufacturing environment supports the stable operation of facilities, smooth management of systems, corrective action and preventive action (CAPA), when deviations occur.

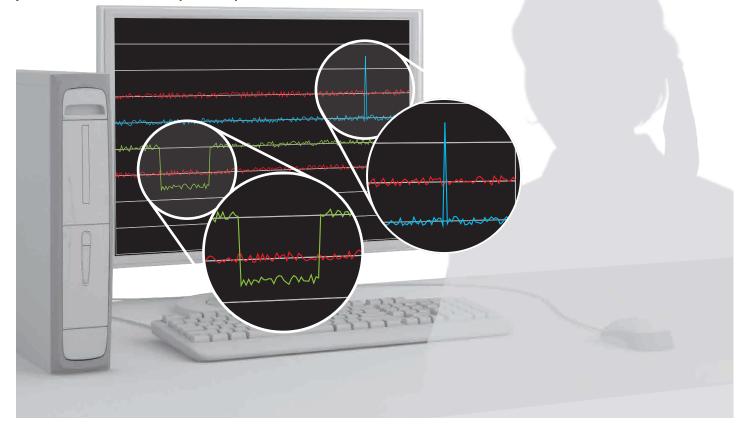



Features of Air-conditioning Monitoring System TAICOM-RP Monitor

■ Continuous Monitoring

The temperature, humidity, room pressure, hydrogen peroxide concentration and etc. can be monitored continuously.

Alarm Monitor


The alarm thresholds can be set individually.

Security

Comply with 21 CFR Part 11. (Electronic signature, alteration prevention, audit trail, etc.) Since operation authority can be added for each user, the loss of data by unexpected operation can be prevented.

CSV

Performs computerized system validation (CSV) based on GAMP5.

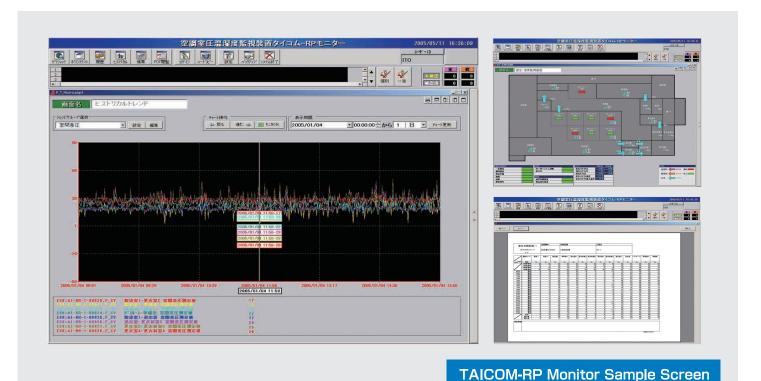

Standard Hardware

① **OS**····· Windows 7 (32bit)

2 Memory 4GB

 \bigcirc HD············· 120GB \sim 1TB

Software

- · iFIX
- · Adobe Acrobat
- · INSTANT COPY
- · SPREAD